


Picarro & Costech Combustion Module – Installation Requirements and Schematic

Items marked in **red** need to be purchased by the customer ahead of an installation. Items marked in **blue** also need to be purchased by the customer, or alternatively they are provided in the optional Picarro O2 and N2 Dry Gas Kit (Picarro P/N: A0924). If additional modifications to the stainless steel tubing are made, please purchase more tubing and additional Swagelok nuts and ferrules (Part: SS-200-NFSET).

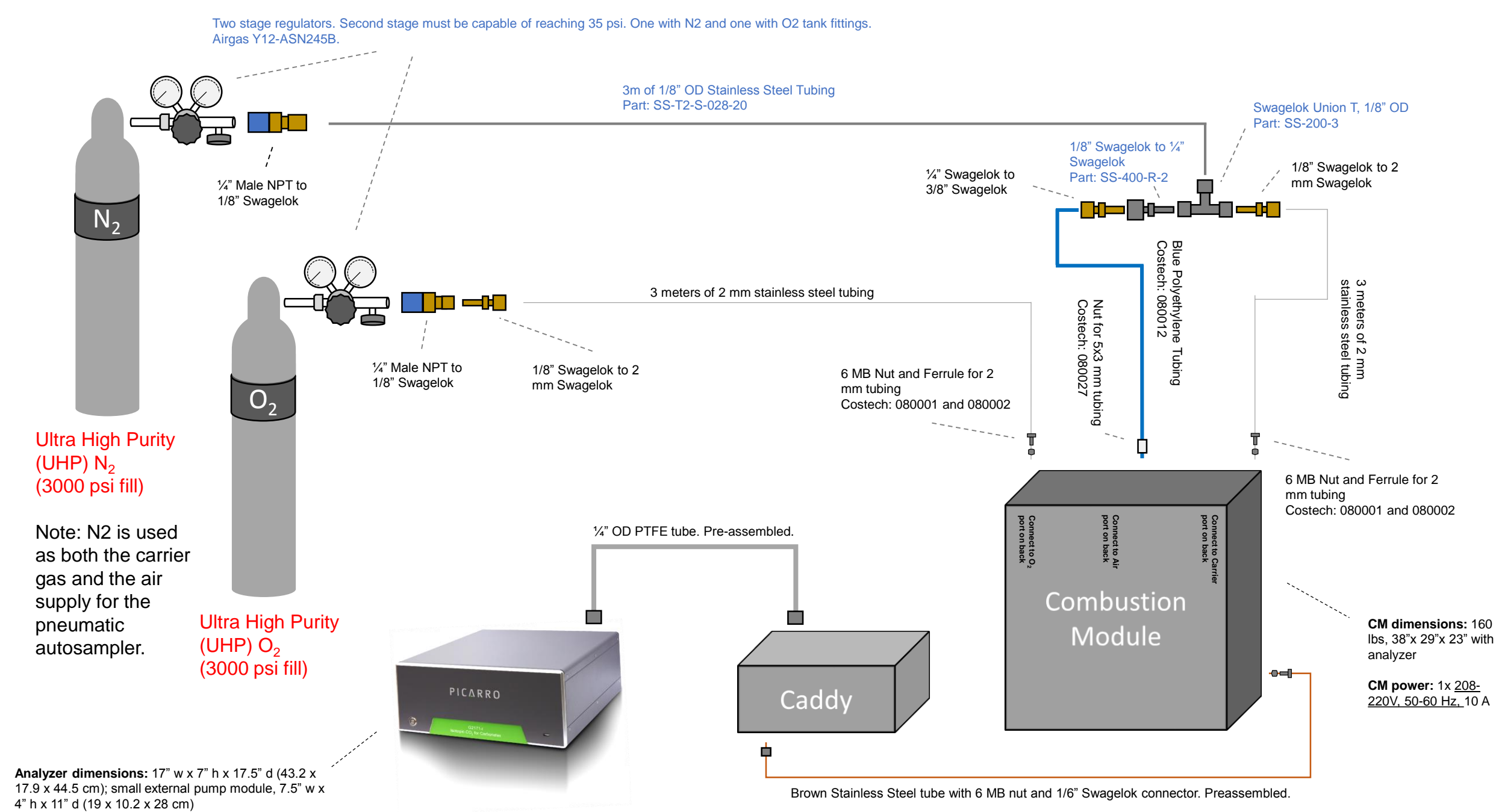

Ultra High Purity (UHP) N₂ (3000 psi fill)

Note: N2 is used as both the carrier gas and the air supply for the pneumatic autosampler.

Ultra High Purity (UHP) O₂ (3000 psi fill)